

**Jazzy®**  
1113ATS

**VERSATILITY IN A COMPACT,  
PORTABLE PACKAGE**

## Jazzy® 1113 ATS

The Jazzy® 1113 ATS offers everything you need—all of the “big chair” options and accessories, even rehab capabilities—on the most compact frame available. You’ll take on tight spaces and narrow hallways without sacrificing a strong outdoor performance never before seen in a chair of this size. For sharp ramp transitions as well as a smooth and stable ride even over exceptionally rough and uneven terrain, the Jazzy 1113 ATS is equipped with Active-Trac® Suspension. As an added plus, the 1113 ATS fully disassembles for transport.

**JAZZY® 1113 ATS**  
shown in Viper Blue with  
medium-back seat, PG VSI  
controller and foot platform



**Pride®**  
Mobility Products Corp.

800-800-8586 (US) • 888-570-1113 (Canada) • [www.pridemobility.com](http://www.pridemobility.com)

JZ 1113 ATS SS 01/07



The Jazzy® 1113ATS series models include:

- Jazzy 1113 ATS 2SP-SS
- Jazzy 1113 ATS 2SP-C

#### FEATURES

- Easily disassembles for transport
- Powerful two-motor, mid-wheel drive design offers outstanding maneuverability
- Low center of gravity design delivers a very smooth and stable ride
- 10" knobby drive wheels and optimum weight distribution provide superior traction
- On-board battery charger
- Fully programmable controller

#### OPTIONS

- TRU-Balance® Power Tilt<sup>1</sup>
- TRU-Balance® Power Recline<sup>1</sup>
- Contoured high-back seat
- Synergy® seating system
- Contoured medium/high-back with solid seat pan
- Cup holder
- Weather cover
- Rear basket
- Oxygen tank holder
- Cane/crutch holder
- Walker holder
- Swing-away joystick mount
- Safety flag
- Off-board charger


1. 200 lbs. weight capacity with TRU-Balance® Systems

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure.

## Jazzy® 1113 ATS

Weight capacity:	300 lbs.
Maximum speed:	Up to 4 mph
Ground clearance:	2.5"
Turning radius:	20.5"
Overall length:	36"
Overall width:	23"
Wheels:	Drive: 10"
Front anti-tips:	6" solid
Rear casters:	8" solid
Suspension:	Active-Trac®
Drive train:	Two-motor, mid-wheel drive
Braking system:	Intelligent braking (electronic, regenerative disc brakes)
Std. Electronics:	50 amp, PG VSI controller
Specialty controls:	Yes
Battery charger:	On-board, 4 amp (standard)
Per-charge range:	Up to 15 miles
Battery requirements:	(2) 12 volt, deep cycle
Battery size:	U-1
Transport weight	24.5 lbs. each
Battery weight:	
Base weight:	
Std. Seat dimensions:	18" x 18"
Warranty:	Lifetime limited warranty on frame; 2-year limited warranty on electronics; 18-months limited warranty on drive motors

Colors:

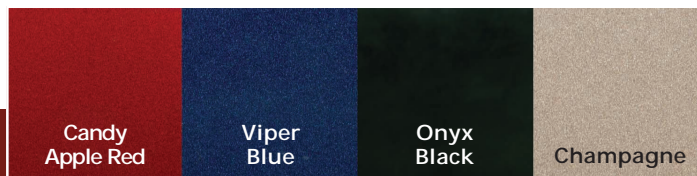

- Full Platform
- Power chair lift with straps
- Backpacker™



Backpacker™



Full Platform

These lifts are compatible with the Jazzy® 1113 ATS.

