

Ti

L I T E®



Aero X Shown with Available Options

Aero X

THE NEW TILITE AERO X. EXTENDING THE
DEFINITION OF THE ULTIMATE RIDE™.

THE NEW TILITE AERO X. EASING EVERYDAY MOBILITY!



The new aluminum Aero X combines unmistakable TiLite design elements with the convenience of a folding chair. Available in two versions, Fixed or Swing In or Out Removable Front Ends, the Aero X is an affordable everyday chair that utilizes quality design and new features to ensure maximum function and durability. By combining TiLite's signature "bullet" caster housing with a unique curved axle plate, the Aero X is as adjustable as you need it to be. Test out the Aero X and you will understand why TiLite is known for building The Ultimate Ride.

TILITE INNOVATION.


- "Tru-Fit" system.
- Rigid Folding Mechanism - "RFM".
- Adjustable "Bullet" caster housing.

SIT DOWN AND STAND OUT. The Aero X features TiLite's new Tru-Fit clamp-on component system; comprised of the RFM, rear axle plate and front caster clamp. These components all clamp onto the new Aero X frame instead of being bolted through. The Tru-Fit system results in a 50% reduction of holes in the frame, allowing for finite adjustments since there is no need to line up components with pre-drilled holes. This design innovation also greatly improves the durability of the Aero X frame.

RFM. TiLite's exclusive Rigid Folding Mechanism RFM cross tube design, first developed for the TX, allows TiLite to incorporate one-piece side frames into its newest line of folding chairs. This innovation increases performance by reducing flex because the side frames are welded together, not bolted. TiLite's RFM cross tubes also feature snap in seat saddles that lock into place, helping the Aero X achieve near rigid chair performance.

COLOR. TiLite in color? TiLite's Aero X is available in 14 frame colors. Whether you are looking for beautiful or just plain cool, now there is a TiLite for you. Mix and match TiLite's frame and accessory colors on the Color Matrix at TiLite.com.

FRAME COLORS

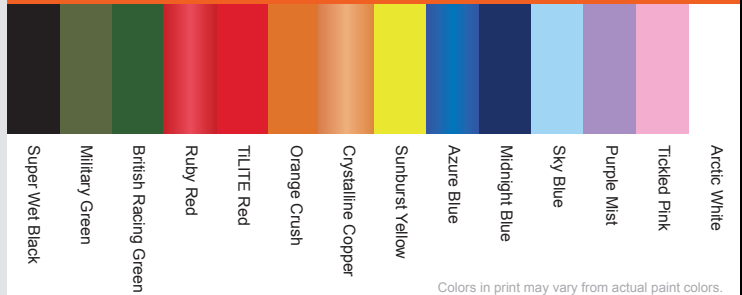

Colors in print may vary from actual paint colors.


TECHNICAL SPECIFICATIONS

Frame Type	Aluminum Folding
Frame Style	Fixed Front End - Full Adjustability Swing In or Out Removable Front End - Full Adjustability
Weight Limit	250 lbs.
Frame Colors	14 Color Options
Seat Width	12" - 20"
Seat Depth	12" - 20"
Front Seat Height	15" - 22"
Rear Seat Height	13" - 20.5"
Camber	0°, 2°, 4°, 6° or 8°
Front Angle	Fixed Front - 70°, 80° or 85° Swing Away - 70°, 80°, 85° or 90°
Front Wheels	18 Options
Rear Wheels	50 Options
Handrims	9 Options
Center of Gravity	1" - 6"
Chair Weight*	Fixed Front End - 16.0 lbs. Swing In or Out Removable Front End - 18.4 lbs.



HAND MADE

PERFORMANCE

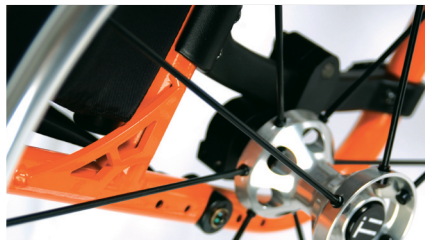

* Weights based on 16" x 16" wheelchairs with all standard features except rear wheels and wheel locks. Rear wheels are not included due to the varying weights of the numerous wheel configurations offered by TiLite. For reference, a pair of standard 24" rear wheels weighs 7.6 lbs. A pair of standard wheel locks weighs 0.9 lbs.



STYLE