



THE ROHO[®] nexus[®] SPIRIT[®] CUSHION

**SETTING THE
NEW STANDARD
IN WHEELCHAIR
CUSHIONING**

The nexus SPIRIT Cushion is unlike any ROHO Cushion you've ever seen before. It combines the stability of a contoured foam base with the therapeutic care for which ROHO[®] DRY FLOATATION[®] cushions are known. The contoured foam base is designed to allow for increased stability for transfers and provides exceptional positioning of the pelvis and lower extremities for enhanced sitting posture. The ROHO Cushion, positioned to address critical pressure points, provides you with the highest level of skin protection from the originators of ROHO DRY FLOATATION Technology.

Life is full of challenges, sitting shouldn't be one of them.



nexus[®] SPIRIT[®] CUSHION

SKIN PROTECTION.

The nexus SPIRIT provides protection to the ischia, sacrum and coccyx area through the placement of a ROHO DRY FLOATATION support pad. Unparalleled in performance, the ROHO Cushion provides an adjustable, low friction/low shear environment that can protect the user against wounds.

STABILITY.

The contoured foam base provides increased stability for transfers in and out of the cushion.

ENHANCED SITTING POSTURE.

The contoured foam is designed to provide improved positioning of the pelvis and lower extremities, enhancing sitting posture.



APPLICATIONS

THE nexus SPIRIT CUSHION IS MOST EFFECTIVE FOR USERS WHO:

- are at moderate risk for skin breakdown
- have chronic pain complaints
- require increased stability for transfers
- have difficulty with lower extremity positioning
- are able to perform independent weight shifts

FEATURES & BENEFITS

LIGHTWEIGHT:

Easy to handle and carry.

WASHABLE:

DRY FLOATATION support pad may be cleaned with soap and water. Foam base may be wiped down. Fluid resistant cover is machine washable.

REPAIRABLE:

Neoprene rubber support pad can be patched or sent to ROHO for any necessary repair.

SPECIFICATIONS

INCLUDES:

nexus SPIRIT cushion, two-way stretch cover, hand inflation pump, repair kit and operating instructions.

CONSTRUCTION:

Flame-resistant neoprene rubber support pad, flame-resistant polyurethane seating foam.

SIZES:

Available to fit chairs: 14" x 14" / 36cm x 36cm, 14" x 16" / 36cm x 41cm, 15" x 15" / 38cm x 38cm, 16" x 16" / 41cm x 41cm, 18" x 16" / 46cm x 41cm, 18" x 18" / 46cm x 46cm and 20" x 18" / 51cm x 46cm.

(NOTE: Sizes listed represent chair size. Actual cushion size is approximately 1.5³/₄ cm deeper to allow for placement between back support canes. Customization only available for neoprene support pad.)

WEIGHT:

3* lbs/ 1 kg (Based on NS1616C model)

*Varies by size

WEIGHT LIMIT:

Unlimited¹

¹Cushion must be properly sized to client.

WARRANTY:

24-Month Limited Warranty

MEDICARE PART B CODE:

K0734 Adjustable Skin Protection and Positioning Wheelchair Cushion

We recommend consulting with a certified rehab professional when ordering and fitting ROHO Group products. Visit NRRTS.org for information on certified providers in your area.

MORE PRODUCTS FROM ROHO




CONTOUR SELECT™ CUSHION



JETSTREAM PRO™ BACK SYSTEM

FOR ADDITIONAL INFORMATION, PLEASE CONTACT

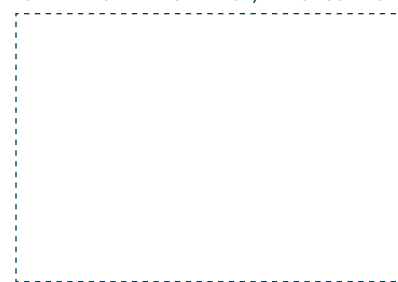

100 North Florida Avenue • Belleville, IL 62221-5429 • USA • www.therohogroup.com
 U.S. & Canada: 1 800 851 3449 • 1 618 277 9173 • fax 1 888 551 3449
 Outside the U.S. & Canada: 1 618 277 9150 • fax 1 618 277 6518

©2006 ROHO, Inc. The following are registered trademarks of The ROHO Group: ROHO®, shape fitting technology®, nexus® SPIRIT®, JetStream Pro™, CONTOUR SELECT™, and DRY FLOATATION®.

ROHO, Inc. has a policy of continual product improvement and reserves the right to amend specifications presented in this sell sheet.

LITG00220
 SCHW 2M 5/08