



Shown with R130 Standing sling

### Features

- Innovative adjustment features adapt to a wide variety of body sizes and shapes. That means comfort at any height, width and weight range up to 350 lb.
- Low friction casters significantly improve rollability and contribute to caregiver safety
- Multi-function slings allow for quick toileting, stand assist or full-support seated transfer
- Adjustable leg support features a soft pad with lateral contour that secures legs in three positions for comfort at any height
- Wide, non-slip footplate provides stable base for transfer
- Under-bed clearance of 4.5" allows accessibility to most nursing home, hospital and home care bed styles
- Available with power opening base (model no. RPS350-2)
- Removable battery eliminates downtime
- Ergonomic features make it ideal for toilet transfers, as well as for dependent and rehabilitation purposes
- Features include fast, simple sling attachment that protects against accidental disengagement; wide push handle; good access for clothing removal; and comfortable slings with extended back support that prevents pressure in the under-arm area
- Adjustable leg strap option secures leg against knee pad for optimum positioning on footplate
- Limited warranty includes three years on lift, one year on motor
- Manufactured in an ISO13485 facility



**Electronics**


- Removable battery box is conveniently located on mast; it features an emergency OFF button and recharge indicator (with audible recharge alarm)
- Quick-release battery eliminates downtime
- Anti-entrapment feature stops motion if boom meets resistance while lowering
- Heavy-duty hand pendant and cord
- Includes two batteries and a recharging station

**Specifications**


Model no.	Pwr. Base RPS350-1	RPS350-2
<b>Height</b>		
Maximum at sling hook-up	63.7"	63.7"
Minimum at sling hook-up	39.6"	39.6"
Base height (clearance)	4.5"	4.5"
<b>Base width</b>		
Open	36.8"	35.8"
Closed	25.7"	25.7"
<b>Base length</b>		
	35.5"	35.5"
<b>Casters</b>		
Front	3"	3"
Rear (locking)	5"	5"
<b>Weight capacity</b>		
	350 lb.	350 lb.
<b>Product weight</b>		
Out-of-carton	115 lb.	118 lb.
Shipping	125 lb.	128 lb.



Model no. RPS350-2 Stand-Up Lift

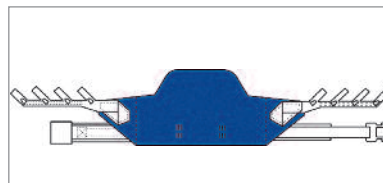

Model no. R130 Extra-large standing sling  
Sling model numbers R130 and R134 are also compatible with the Reliant 440 lift. Complete instructions for lift usage and sling procedures are found in the Invacare operating manual. Only Invacare® slings and accessories should be used on Invacare® lifts.

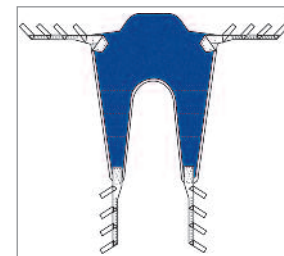

Model no. R134 Extra-large transfer sling

**Sling Specifications**

	R130	R134
Width	60"	34"
Overall Length	16"	42"
<b>Product weight</b>		
Out-of-carton	1.1 lb.	1.3 lb.
Shipping	1.6 lb.	1.8 lb.