

Advanced Pressure Redistribution System



THE ANSWER TO THE CHALLENGES OF COST-EFFECTIVE THERAPY

In today's challenging health care arena, never before has the need for a clinically proven support surface that can help prevent and manage skin breakdown been so important.

Who better to provide you with a system you can trust, than the company that has been the leader in cost-effective support surfaces for more than 40 years?

The ProActiveAir Pressure Redistribution System is a proven winner and can provide your residents with incredible care and protect your facility in these very litigious times.




Eight anatomically inspired and alternatively connected air sectors provide both clinical support and “exchange” of air as needed and provides a custom pressure redistribution profile for residents of almost every body size, shape and weight.

Unique, multi-zoned topper layer that provides incredible comfort.

Airport for easy connection to a control unit to provide vital alternating pressure therapy.

Heel slope for added protection for the critical heel area with contour cuts for comfort — allows the “meatier” calf to bear more of the weight of the lower leg.

Pressure Mapping 6'1"/210 lb 49-year-old male subject

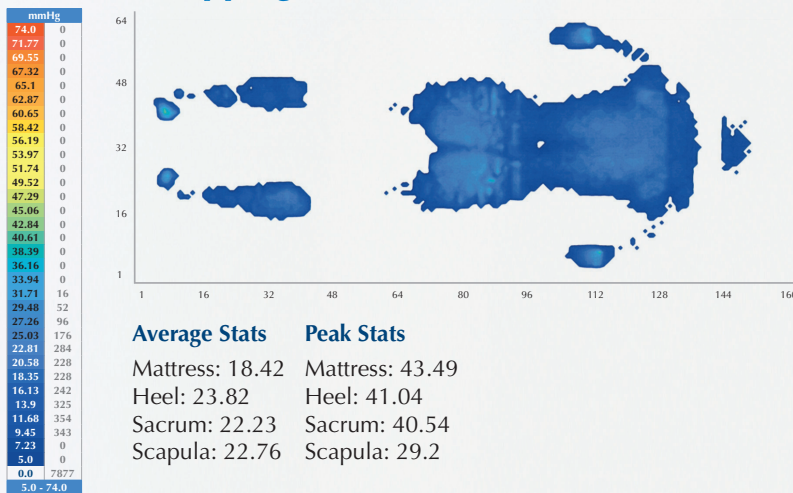

System Features

- Four-way stretch top cover — waterproof with very low MVTR
- Flexskin inner cover (not shown in photograph) — provides a buffer under the top cover to help lower friction and shear forces
- Firm perimeter foam borders for stability during ingress/egress — includes articulation channels for bed articulation
- Evacuation handles
- 360° zipper for easy access to the inner core
- FlameTex Fire Protection to meet all applicable fire codes including 16CFR1633 (not shown in image)

Details

- Weight capacity: 500 lbs
- Warranty
 - Mattress: 5-year limited
 - Control Unit: 1-year (sold separately)
- Mattress components are treated with an anti-bacterial/anti-microbial to help prevent bacterial growth and resident infections
- Available in normal bed dimensions (35 x 75, 80, 84). Custom sizes are available at additional cost

Available Options

- Resident Reminder System (RRS) — raised side borders to help prevent accidental falls
- A variety of top and bottom covers and colors including a non-slip bottom cover

Ordering Information

Product	Size	Order Number
The ProActiveAir Advanced Pressure Redistribution System (with Airport for Control Unit Connection)	35 x 75 35 x 80 35 x 84	AT-1BAP1J-675 AT-1BAP1J-680 AT-1BAP1J-684
Alternating Pressure Control Unit	35 x 84	ECU1
Premium Alternating Pressure Control Unit	35 x 75	ECU2