



## NEW! THE ROHO® HYBRID ELITE® CUSHION

### POSITIONING

### STABILITY

### SKIN PROTECTION

ROHO® offers you more choices. Introducing the new ROHO Hybrid Elite® Cushion – the wheelchair cushion that combines the superb skin protection and shape matching ability of ROHO DRY FLOATATION® with the stability of Jay® base technology. This new cushion design allows for improved pelvic alignment and stability for improved position and function while maintaining optimal immersion and envelopment of both the ischial tuberosities and greater trochanters. The ROHO shape fitting cushion insert provides critical pressure redistribution around the pelvis area providing you with the superior skin protection that ROHO is known for. Simply adjust the cushion insert to fit the client's shape in order to minimize tissue deformation and maximize blood flow.

Skin protection, positioning and stability.

**Two trusted brands. One effective cushion.**



# NEW! ROHO® HYBRID ELITE® CUSHION

## POSITIONING.

Achieve greater ability to position your client with the ROHO Hybrid Elite cushion. Its customized Jay contoured foam base provides you with the opportunity for both lower extremity and pelvic alignment and control.

## SKIN PROTECTION.


The ROHO Hybrid Elite cushion hydrostatically loads the ischial tuberosities and the posterior shelf of the greater trochanters for those clients who are at high risk of breakdown. Unparalleled in performance, the ROHO cushion insert allows for immersion and envelopment of the client's individual shape. It provides an adjustable, reduced friction and shear environment that can help protect the client from acquiring pressure ulcers from sitting. The dual valve option on the ROHO insert allows for additional customization to match a client's specific shape and asymmetry.

## STABILITY.

The ROHO Hybrid Elite's customized Jay contoured foam base provides improved positioning of the pelvis and lower extremities, enhancing sitting posture. This added stability assists the client in achieving a safer transfer in and out of their wheelchair.



ROHO combines a customized Jay contoured foam base for increased stability and a ROHO DRY FLOATATION Cushion Insert for skin protection. Available in single and dual compartment models.

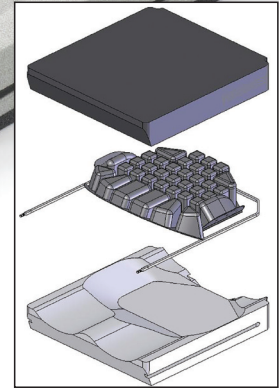

## APPLICATIONS

### THE ROHO HYBRID ELITE CUSHION IS RECOMMENDED FOR CLIENTS WHO:

- are at a moderate to high risk for skin breakdown
- require pelvic and lower extremity positioning
- have difficulty with lower extremity positioning and alignment
- have asymmetrical postures and need front to back or side to side positioning (dual compartment sizes applicable)

## FEATURES & BENEFITS

### ADJUSTABLE:

Cushion insert matches the individual's shape by adjusting the amount of air.

### WASHABLE:

Air cushion insert is easily cleaned with soap and water. Cover with encased comfort foam is machine washable. Foam base may be wiped with damp cloth.

### REPAIRABLE:

Neoprene rubber cushion insert can be patched or repaired to extend use.

### FLUID RESISTANT:

Customized Jay foam base and ROHO air insert are fluid resistant.

We recommend consulting with a certified rehab professional when ordering and fitting ROHO products. Visit [NRRTS.org](http://NRRTS.org) for information on certified providers in your area.

## SPECIFICATIONS

### INCLUDES:

ROHO cushion insert, customized Jay contoured foam base, two-way stretch cover, hand inflation pump, repair kit and operation instructions

### CONSTRUCTION:

ROHO cushion insert: flame resistant, neoprene rubber  
Customized Jay contoured foam base: flame resistant, fluid resistant closed cell foam

### SIZES:

Wide Range, to fit chair widths 14" to 24"/35cm to 60cm

### WEIGHT:

Varies by size; approximately 3.5 lbs. / 1.5 kg (based on 17.75 x 18.5 size)

### WARRANTY:

24-month Limited Warranty

### WEIGHT LIMIT:

500 lbs / 225 kg

### MEDICARE PART B CODE:

K0734/5 Adjustable Skin Protection Wheelchair Cushion

## MORE PRODUCTS FROM ROHO




### HIGH PROFILE® QUADRO SELECT®



### JETSTREAM PRO® BACK SUPPORT SYSTEM

### FOR ADDITIONAL INFORMATION, PLEASE CONTACT:

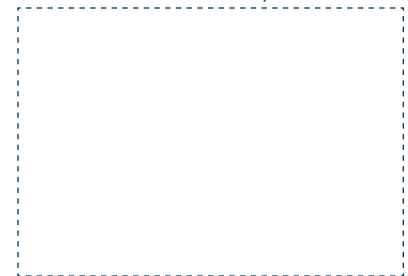

100 North Florida Avenue • Belleville, IL 62221-5429 • USA • [www.therohogroup.com/hybrid-elite](http://www.therohogroup.com/hybrid-elite)  
U.S.: 1 800 851 3449 • 1 618 277 9173 • fax 1 888 551 3449  
Outside the U.S.: 1 618 277 9150 • fax 1 618 277 6518

©2009 ROHO, Inc. The following are trademarks and registered trademarks of ROHO, Inc.: ROHO®, shape fitting technology®, Hybrid Elite®, HIGH PROFILE®, JET STREAM PRO®, DRY FLOATATION®, and QUADRO SELECT®. The following is a registered trademark of Sunrise Medical, Inc.: Jay®.

ROHO, Inc. has a policy of continual product improvement and reserves the right to amend specifications presented in this sell sheet.

LITG02132  
Rev. 09/10 10/10